



Volumetric Flasks

Class A, Low-Actinic Amber, F Stopper, To Contain

SIBATA 2307A-Series Low-Actinic Amber Volumetric Flasks meet ASTM E-288 Class A Specifications & USP Standards for Volumetric Glassware.

These flasks feature machine-blown bodies and heavy-beaded, heavy-tubing necks, tooled for F stoppers. The graduation line is sharp and permanent; large, white, block-letters allow the labeling to be easily read. Each flask also has a large marking square. The low-actinic amber stained glass provides protection for materials sensitive to light.


All sizes are supplied with a glass F stopper.

Specifications

Cat. No.	Capacity (mL)	Tolerance (\pm mL)	F Stopper No.	Approx. O.D. x Height (mm)
2307A-10	10	0.02	9	28 x 95
2307A-20	20	0.03	9	36 x 115
2307A-25	25	0.03	9	39 x 122
2307A-50	50	0.05	13	49 x 150
2307A-100	100	0.08	13	60 x 180
2307A-200	200	0.10	16	75 x 220

Packaging

Cat. No.	Pack Qty.	Case Qty.
2307A-10	6	72
2307A-20	6	72
2307A-25	6	72
2307A-50	4	48
2307A-100	4	48
2307A-200	4	24

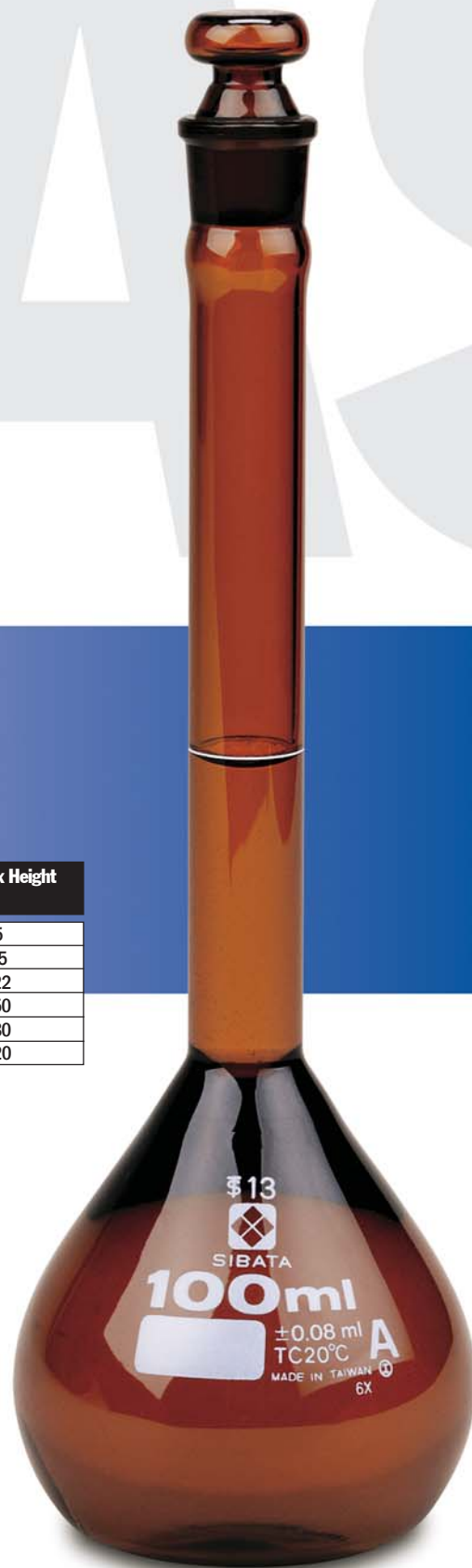

VEE GEE Scientific, Inc.

13600 NE 126th Place, Suite A • Kirkland, WA 98034
(800) 423-8842 • (425) 823-4518 • Fax: (425) 820-9826
www.veegee.com • e-mail: sales@veegee.net